



Hillrom™


Envision® Therapy Surface
Rental offering

ADVANCED WOUND THERAPY

The Envision® Therapy surface represents one of our most advanced and complete powered air surface solutions, promoting better outcomes and enhanced comfort for patients with high pressure injury risk.

SUPPORTING SAFETY AND CAREGIVER EFFICIENCY

- Head-of-bed angle indicators and alert
- Turn Assist
- Bed Exit
- Patient weight capacity: 400 lbs
- Touchscreen Display

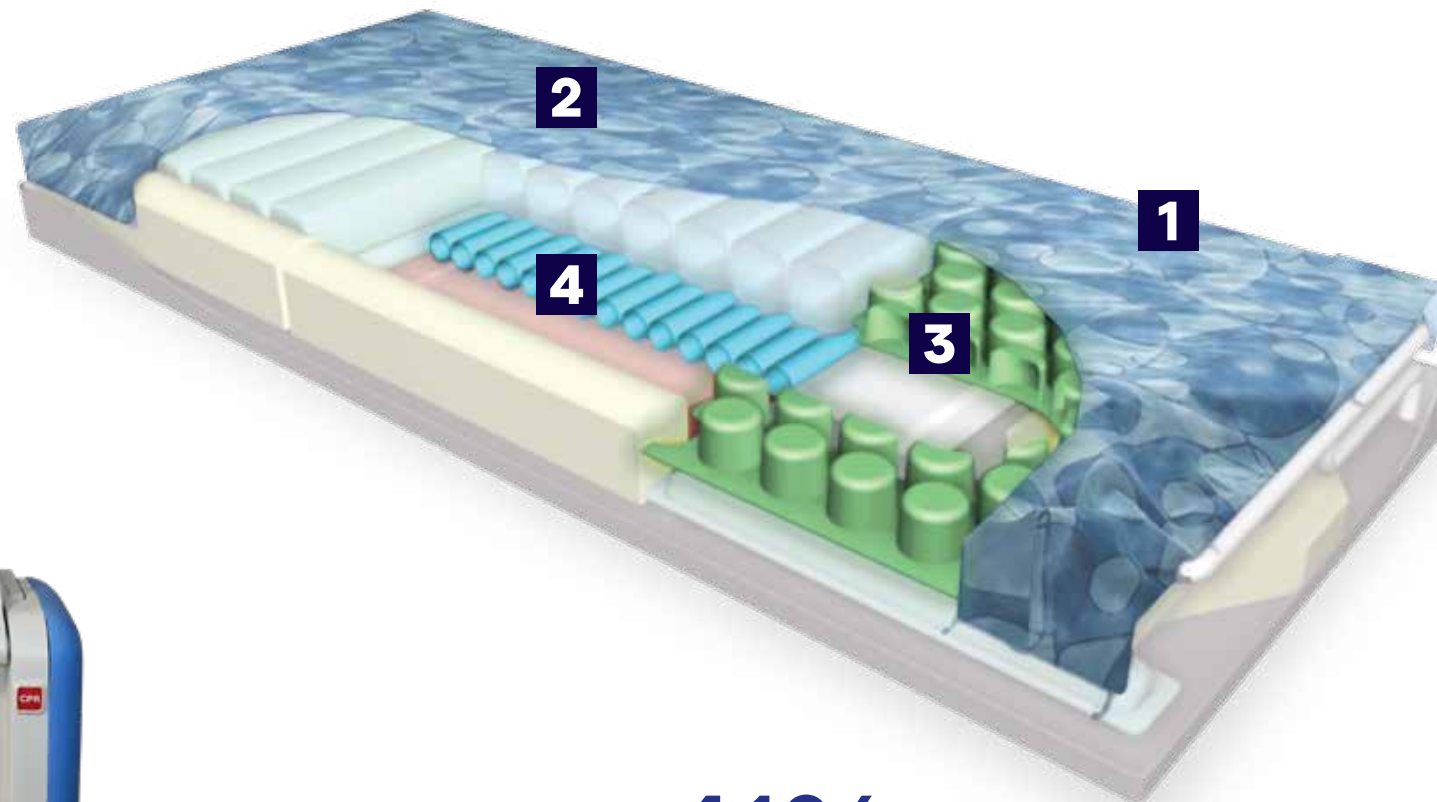

41%↓


DECREASE IN PRESSURE INJURY VOLUME

in seven days or less, according to a 30-patient ICU step-down study.¹

- 1 ADVANCED MICROCLIMATE® TECHNOLOGY**
Technology designed to push airflow underneath the patient from one end of the surface to the other, encouraging more removal of heat and moisture.
- 2 SHEARRELIEF® SYSTEM**
The ShearRelief® Algorithm feature automatically adjusts cushion pressures when the head-of-bed is raised. This, along with non-crushable 3-D fabric that shifts with the patient, helps minimize shear and friction.
- 3 VERTICAL AIR CELL TECHNOLOGY**
Multiple vertical air cells act independently, optimizing immersion and envelopment of the vulnerable heel area.
- 4 PRESSURESMART® SENSING TECHNOLOGY**
A highly responsive sensor system reacts to the patient's position in real-time, adjusting air cushions to minimize pressure.

Advanced Microclimate® Technology

Advanced Microclimate technology pushes airflow underneath the patient from one end of the surface to the other. This feature, combined with breathable, crush-resistant fabric in the top liner, encourages more thorough removal of heat and moisture.

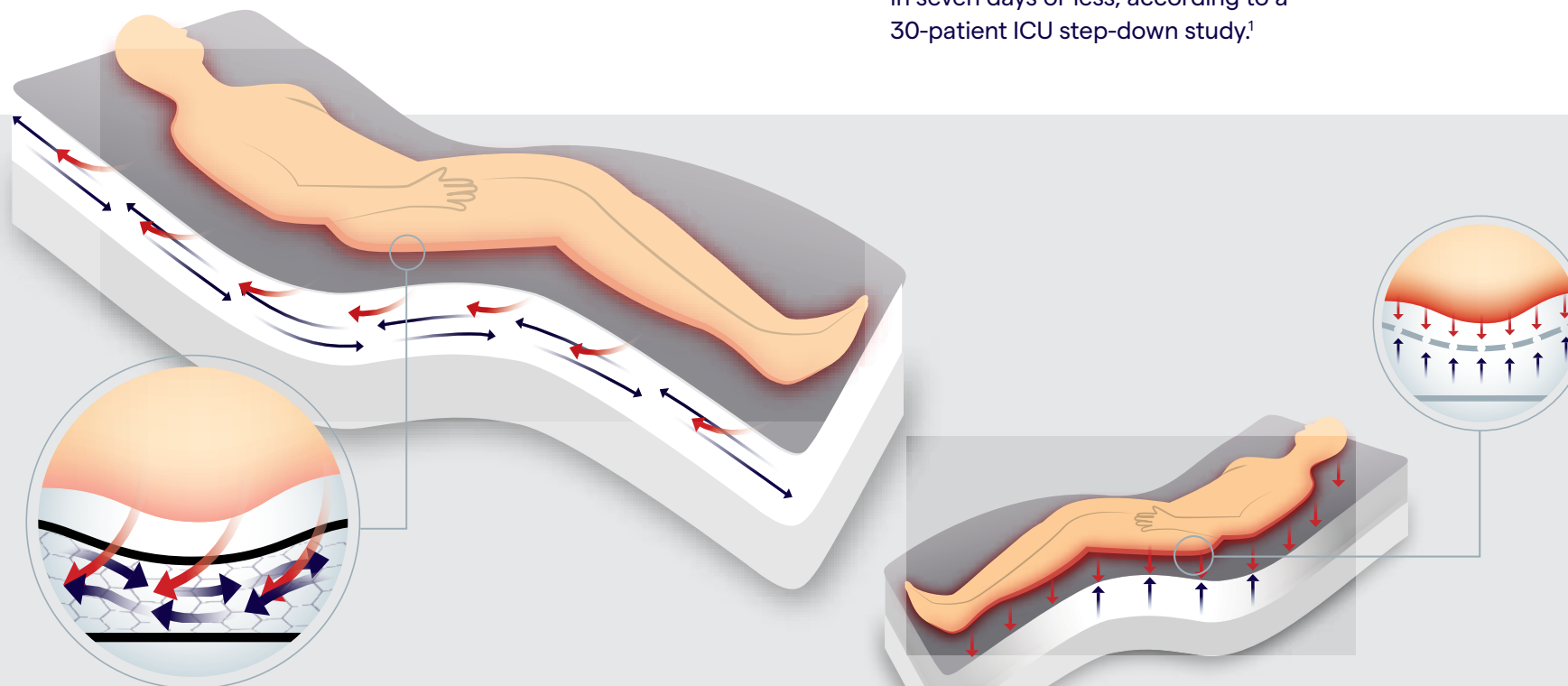

TRADITIONAL LOW AIR LOSS

Traditional low air loss systems have pinholes to provide airflow, potentially leaving more heat and moisture beneath the patient.

TECHNICAL SPECIFICATIONS

Mattress Width (inflated)	35.5" ± 0.5" (90.2 cm ± 12.7 mm)
Mattress Length (inflated)	84.5" ± 0.5" (214.6 cm ± 12.7 mm)
Mattress Height (inflated at maximum pressures for a step deck bed such as the VersaCare® Bed or TotalCare® Bed System)	10.875" ± 1" (27.623 cm ± 3 cm)
Mattress Height (inflated at maximum pressures for a flat deck bed such as the Advanta™ 2 Bed or CareAssist™ Bed)	8" ± 1" (20 cm ± 3 cm)
Air Supply Unit Dimensions	14.5" W x 13.75" H x 8.5" D
Power Cord Length	168" +/- 1" (427 cm +/- 3 cm) See user manual
Power Requirement	120 V AC, 60 Hz, 3A
Recommended Therapeutic Patient Weight Range	70 lbs to 400 lbs (32 kg to 181 kg)



Multi-Frame Compatibility

As a therapy mattress replacement system, the Envision surface is versatile. It's compatible with many Hillrom™ frames, including the VersaCare, TotalCare, Advanta 2 and CareAssist beds.



Hillrom™

For more information, please contact your local Hillrom sales representative at 1-800-445-3730 or order online at rental.hillrom.com.

hillrom.com

¹ Fimiani J, Cohen E, VanGilder C. A pilot study of a new low air loss treatment surface in the critical care setting. Poster Presentation. Symposium on the Advances of Skin and Wound Care; April 25-29, 2008; San Diego, CA.

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.