



**Standard Features**

**ADVANTA™ 2 BED**

<b>Emergency CPR:</b> Release lever with dampened descent	+
<b>Brake Alarms:</b> Alerts caregiver when the bed is plugged in and the brake is not set	+
<b>Chair Position:</b> Caregivers can place patient in a chair position with the touch of one button; bed remains at a low height	+
<b>Battery Backup:</b> Backup for all bed articulation (high/low position, head-of-bed, knee section and foot angle)	+
<b>Point-of-Care® Siderail Controls:</b> All caregiver controls are conveniently located in a central location, on both sides of the bed, where patient and caregiver interaction occurs	+
<b>Line-of-Site® Feature:</b> Head-of-bed and Trendelenburg angle indicators integrated into siderails	+
<b>Electric Vascular Foot:</b> Allows the caregiver to raise patient’s feet without having to lift the weight of the bed or patient’s legs	+
<b>OneStep® Siderail Release Mechanism:</b> One-handed up/down siderail operation	+
<b>Built-in Bed Extender:</b> Easily adds length to the bed	+
<b>Open Architecture (Range of available mattresses)</b>	+
<b>SlideGuard® Technology:</b> SlideGuard® bed articulation technology is designed to support patient positioning	+


**Optional Features**

<b>IntelliDrive® Powered Transport System:</b> Allows for easy transport of a patient in bed	<b>Optional</b>
<b>Three-Mode Bed Exit System (Patient Position Monitor):</b> Bed exit system with three different sensitivities: Patient Movement Mode, Bed Exiting Mode, and Out-of-Bed Mode	<b>Optional</b>
<b>Ability to send data to Nurse Call, EMR or other Clinical Applications</b>	<b>Optional</b>
<b>SafeView® System:</b> Safety alerts at the foot end of the bed display green or amber light to show the bed status (brakes not set, bed exit, bed not in lowest position, siderails down); lights flash an amber color alert when bed is not set in any of the programmed statuses	<b>Optional</b>
<b>Integrated Scale</b>	<b>Optional</b>

Hill-Rom®  
**Advanta™ 2 Bed**  
**Spec Sheet**



**Advanta™ 2 Bed – Available Wood Grain and Graphic Finishes**

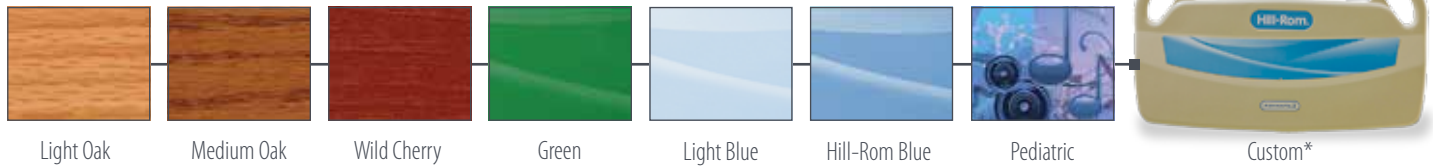

\*Custom options available. Contact your Hill-Rom representative for details.

Technical Specifications		Optional scale system	
Frame width	40" (101.6 cm)	Accuracy	+/- 2.2 lbs or 1.0% of patient weight whichever is greater
<b>Frame length with head/footboard (w/ bumpers)</b>		Bed can weigh in any articulated deck position	Yes
Normal	91" (231 cm)	Scale repeatability	0.2% of patient weight
Extended	94" (238.9 cm)	Low position floor to deck	15.75" (42.55 cm)
<b>Sleep deck</b>		High position floor to deck	32.5" (82.55 cm)
Width	36" (91.44 cm)	Max head elevation	65°
Normal length	84" (213.4 cm)	Max Trend and Reverse Trend	16°
Extended length	87" (221 cm)	Caster size	6" (15 cm)
Bed clearance under frame	5.5" (13.97 cm)	Safe working load	500 lbs

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.

©2013 Hill-Rom Services, Inc. ALL RIGHTS RESERVED.

ORDER NUMBER 183637 rev 2 E 01/07/14 ENG – US

Zinc Number US-AC-0713-0333(1)

[www.hill-rom.com](http://www.hill-rom.com)

USA 800-445-3730  
 CANADA 800-267-2337

Enhancing outcomes for patients and their caregivers:

