



The Volara™ System Oscillation & Lung Expansion Therapy


MAXIMIZE OUTCOMES

The Volara System delivers trusted OLE therapy, combining Continuous Positive Expiratory Pressure (CPEP), Continuous High Frequency Oscillation (CHFO) and nebulizer in one portable, lightweight device.

- Innovative digital sensing technology eliminates the need for Venturi, while achieving exact pressure setting and regulation
- Portable, lightweight, internal digitally-controlled air source enables use throughout the hospital
- Touchscreen interface allows precise setting of pressures and programs
- One-touch access to frequently used presets supports standardization of therapy across patient groups



- Reinforced pole design
- Attached wire storage basket
- Lower base improves stability
- Locking casters for safety

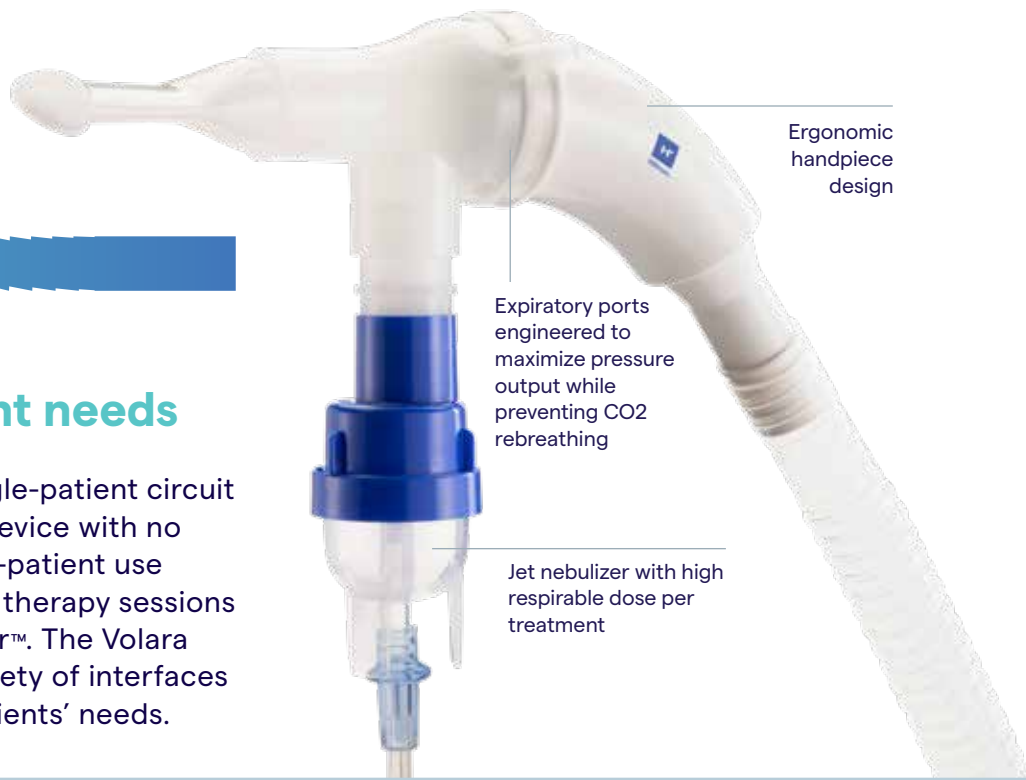

Ergonomic handpiece design


Expiratory ports engineered to maximize pressure output while preventing CO2 rebreathing

Jet nebulizer with high respirable dose per treatment

Meet diverse patient needs

The Volara™ System uses a single-patient circuit that connects directly to the device with no calibration needed. This single-patient use circuit is intended for up to 90 therapy sessions and tracked by our Smart-Filter™. The Volara System can be used with a variety of interfaces and adapters to meet your patients' needs.

- 1** Handset
- 2** Breathing hose
- 3** Nebulizer
- 4** Nebulizer tubing
- 5** Clear spontaneous breathing adapter
- 6** Blue ventilator adapter
- 7** In-line ventilator adapter (22mm x 15mm)
- 8** Mouthpiece
- 9** Tracheostomy adapter
- 10** Adapter (22mm x 22mm)
- 11** O2 Bleed-in adapter
- 12** Smart-Filter with integrated RFID tag

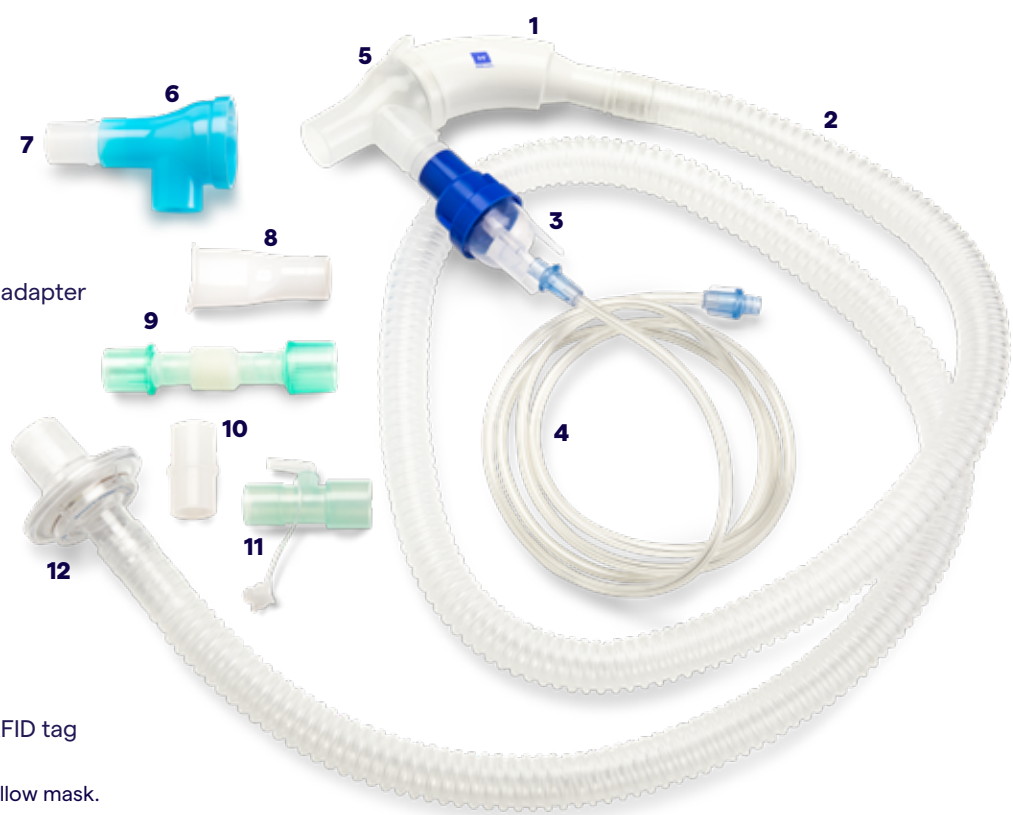

Volara is also indicated for use with a pillow mask. Mask is not included in the circuit kit.

For further information about products or services, please contact your Hillrom sales representative or Hillrom customer service at 1-800-426-4224.

respiratorycare.hill-rom.com

Volara™ and Smart-Filter™ are trademarks of Hill-Rom Services PTE Ltd. Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.