



## The Volara™ System Oscillation & Lung Expansion Therapy


# MAXIMIZE OUTCOMES

The Volara System delivers trusted OLE therapy, combining Continuous Positive Expiratory Pressure (CPEP), Continuous High Frequency Oscillation (CHFO) and nebulizer in one portable, lightweight device.

- Innovative digital sensing technology eliminates the need for Venturi, while achieving exact pressure setting and regulation
- Portable, lightweight, internal digitally-controlled air source enables use throughout the hospital
- Touchscreen interface allows precise setting of pressures and programs
- One-touch access to frequently used presets supports standardization of therapy across patient groups



- Reinforced pole design
- Attached wire storage basket
- Lower base improves stability
- Locking casters for safety

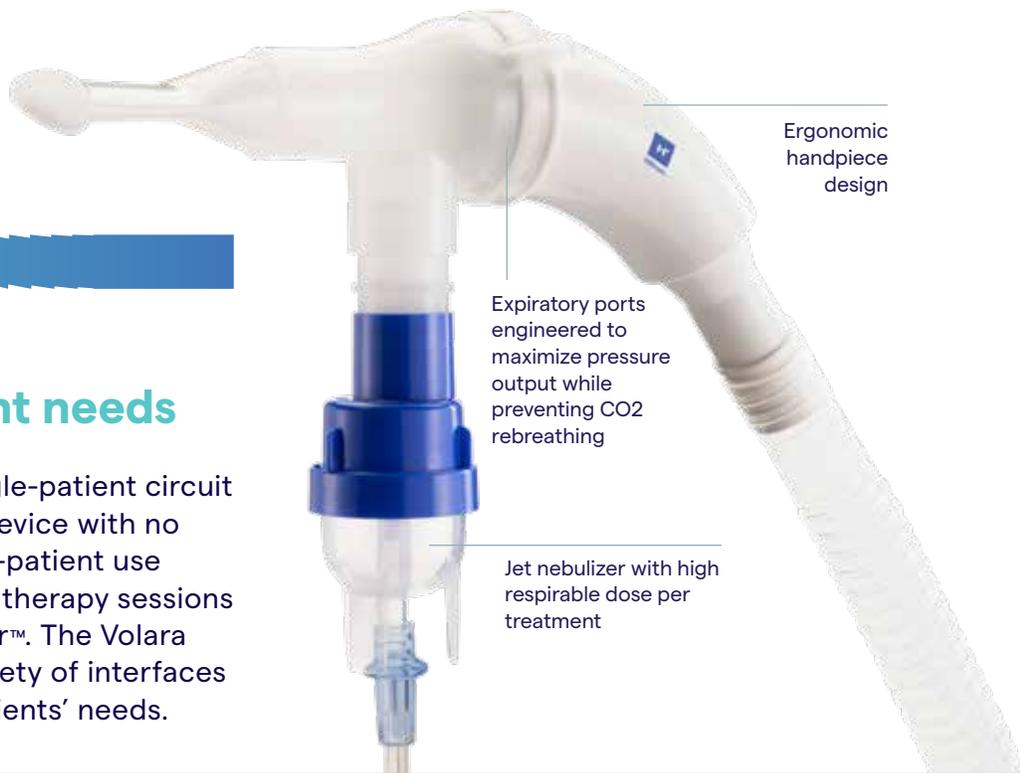


Ergonomic handpiece design

Expiratory ports engineered to maximize pressure output while preventing CO2 rebreathing

Jet nebulizer with high respirable dose per treatment

## Meet diverse patient needs

The Volara™ System uses a single-patient circuit that connects directly to the device with no calibration needed. This single-patient use circuit is intended for up to 90 therapy sessions and tracked by our Smart-Filter™. The Volara System can be used with a variety of interfaces and adapters to meet your patients' needs.

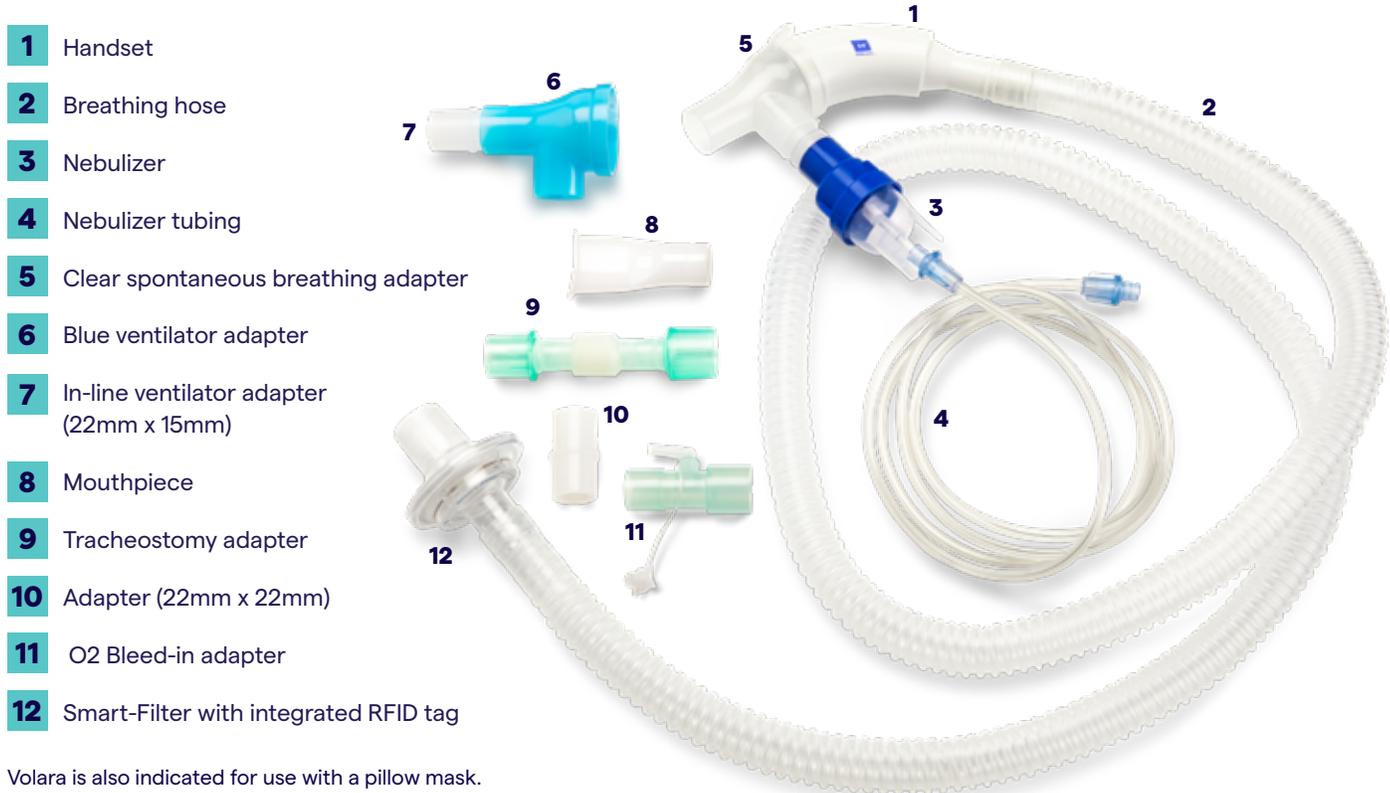

- 1** Handset
- 2** Breathing hose
- 3** Nebulizer
- 4** Nebulizer tubing
- 5** Clear spontaneous breathing adapter
- 6** Blue ventilator adapter
- 7** In-line ventilator adapter (22mm x 15mm)
- 8** Mouthpiece
- 9** Tracheostomy adapter
- 10** Adapter (22mm x 22mm)
- 11** O2 Bleed-in adapter
- 12** Smart-Filter with integrated RFID tag

Volara is also indicated for use with a pillow mask. Mask is not included in the circuit kit.

For further information about products or services, please contact your Hillrom sales representative or Hillrom customer service at 1-800-426-4224.

[respiratorycare.hill-rom.com](http://respiratorycare.hill-rom.com)

Volara™ and Smart-Filter™ are trademarks of Hill-Rom Services PTE Ltd. Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.