




**Hillrom™**

**TotalCare® Bariatric Plus Bed**

# INNOVATIVE SUPPORT TO HELP MOBILIZE PATIENTS OF SIZE WITHOUT COMPROMISING DIGNITY

Many bariatric patients can lose their ability to mobilize or perform self-care activities in only a few days if they do not retain mobility, contributing to complications during their hospitalization.<sup>1</sup>

Use of protocols with early and progressive mobilization of bariatric patients following surgery has been shown to be safe and effective resulting in acceptable LOS and low re-admission rates.<sup>2</sup>

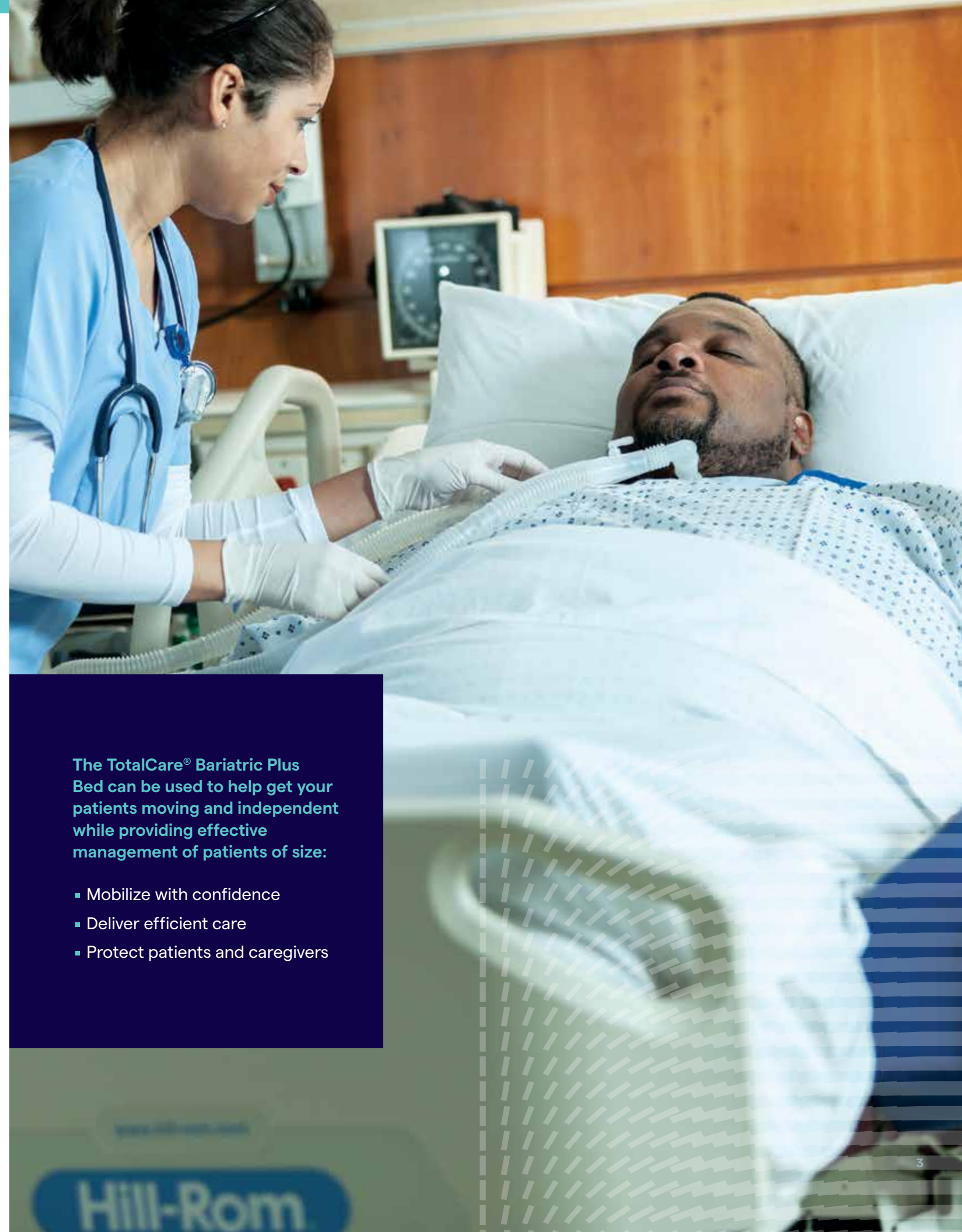

The TotalCare® Bariatric Plus Bed can be used to help get your patients moving and independent while providing effective management of patients of size:

- Mobilize with confidence
- Deliver efficient care
- Protect patients and caregivers

Hill-Rom



## MOBILIZE WITH CONFIDENCE

The features of the TotalCare® Bariatric Plus Bed, together with a Progressive Mobility® Protocol, can provide you with the support you need at the most critical moments to confidently and safely progress your patients from lying in bed, to sitting upright, to standing up.

## DELIVER EFFICIENT CARE

The TotalCare® Bariatric Plus Bed is designed to help you efficiently deliver routine care while comfortably supporting your patients.



- IntelliDrive® powered transport helps enable safe and efficient transport of patients throughout the hospital.
- Turn Assist helps make it easier to perform linen changes and wound checks.
- FlexAfoot™ length adjustment enables quick and easy adjustment of bed to fit patient, minimizing migration and need for repositioning.

# PROTECT PATIENTS AND CAREGIVERS

The TotalCare® Bariatric Plus Bed is designed to support patient safety protocols for helping prevent hospital-associated complications, such as pressure injuries, patient falls and pulmonary conditions related to immobility.

## PRESSURE ULCER MANAGEMENT

- 40" pressure redistribution air surface automatically adjusts to patient's weight, body type, movement and bed position to help treat and prevent pressure injuries.
- Advanced Microclimate® technology helps keep patient's skin cool and dry to protect against skin breakdown.

## FALL PREVENTION

- Bed exit helps supports fall prevention protocols.
- NaviCare® Nurse Call and clinical status board can help identify fall risks.

## PULMONARY MANAGEMENT

- Continuous Lateral Rotation Therapy (CLRT) and Percussion and Vibration Therapy support protocols for helping prevent ventilator-associated complications.
- Head-of-bed angle indicator and alert provide responsive monitoring and protocol compliance.



## TECHNICAL SPECIFICATIONS

<b>Bed</b>	
Low Position - Sleep Deck to Floor	19" (48.2 cm)
High Position - Sleep Deck to Floor	38" (94 cm)
Maximum Overall Width	44.5" (113 cm)
Maximum Overall Length	81.5"-93.5" (206.8-237.3 cm)
Patient Weight Capacity	500 lbs (227 kgs)
<b>Bed Angles</b>	
Head Section, Knee Section	0°-65°, 0°-10°
Max Trend/Reverse Trend	15°/20°
<b>Chair Position</b>	
Head, Thigh, Foot	50°-65°, 0°-10°, 30°-70°
<b>Egress Position</b>	
Head, Thigh, Foot	65°-75°, 0°-10°, 70°-85°
<b>Sleep Surfaces</b>	
Pressure redistribution foam	72"-84"L x 38"W x 7"H (182.7-213.4 cm x 96.5 x 17.8 cm)
Low Air Loss	72"-84"L x 40"W x 11"H (182.7-213.4 cm x 102 x 27.9 cm)
Low Air Loss with CLRT and P&V Modules	72"-84"L x 40"W x 11"H (182.7-213.4 cm x 102 x 27.9 cm)
<b>Caster &amp; Braking System</b>	
Caster Diameter	5" (12.7 cm)
Four Wheel Braking	Yes
Steer and Lock Both Rolling and Swivelling	Yes
<b>Scale System</b>	
Accuracy	1% of patient weight (up to 400 lbs)
Scale Weight Capacity	500 lbs (227 kgs)
<b>Safety Listing</b>	UL60601-1
Power Requirements	120v, 60 cycle AC
Fuse, Maximum Current	15A, 9.9A
<b>Fire Safety</b>	
Meets the requirements of Foam and air surface meet	16CFR1632 CAL129, CAL603, NFPA101, BFDIX-11

## STANDARD FEATURES

- Point-of-Care® Siderail Controls
- Color Touch Screen Graphic Caregiver Interface (GCI)®
- Head-of-bed Alarms (30°/45°)
- Head and Trend Angle Indicators
- Bed Exit
- Battery Backup
- In-bed Scale
- Radiolucent Sleep Deck
- FlexAfoot™ Mechanism

## OPTIONAL FEATURES

- IntelliDrive® Powered Transport (standard on rental)
- Patient Pendant
- Smart Bed Ready – SideCom® System Ready (standard on rental) Nurse Call, Universal TV/Radio, and Lighting (capital purchase only)
- Patient-Helper Trapeze (standard on rental)
- Low Air Loss Surface with Turn Assist
- Low Air Loss Surface with CLRT and Percussion & Vibration Modules



**Hillrom™**

For more information, please contact your local Hillrom sales representative at 1-800-445-3730 or order online at [rental.hillrom.com](http://rental.hillrom.com).

[hillrom.com](http://hillrom.com)

<sup>1</sup> Hahler, Safe Patient Handling and Movement, A Guide for Nurses and Other Health care Providers. 2005. pg 14.

<sup>2</sup> Awad S. Enhanced Recovery After Bariatric Surgery (ERABS): Clinical Outcomes from a Tertiary Referral Bariatric Centre. OBES SURG (2014) 24:753-758.

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.